Therabody

# RecoveryAir Prime

Quick Start Guide

Simplified Recovery.

## Table of Contents

EN	2-4
DE	5-7
FR	8-10
IT	11-13
ES	14-16
NL	17-19
PL	20-22
EL	23-25
NO	26-28
SV	29-31
PT	32-34
RU	35-37
AR	38-40
TR	41-43
SR	44-46
TC	47-49
SC	50-52
KR	53-55
JP	56-58

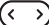
# RecoveryAir Prime

Advanced compression simplified for everybody.

Your RecoveryAir Prime pneumatic compression system comes with a pump, a set of compression boots, blocker plugs, a DC power adapter, and a carrying pouch.


## Feature Callouts
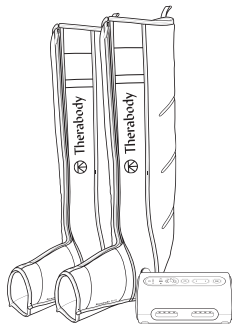
### Buttons


- A**  **Power Button**
- B**  **Start/Pause/Stop Button**
- C**  **Time Button:**  
4 time intervals  
(20, 40, 60, and continuous)
- D**  **Pressure UP/Down Button:**  
  
Adjustable Pressure Range  
(20-100 mmHg) in increments  
of 10 mmHg  
  
Adjustable pressure changes  
via the Therabody app in  
increments of 5 mmHg

### E Indicators

-  **Warning Indicator**  
Indicator lights up solid ON or blinking when  
a device malfunction is detected
-  **Bluetooth Indicator**  
Indicator turns ON only when connected to  
the App
-  **Battery Indicator**  
3 LED lights to indicate battery life

EN

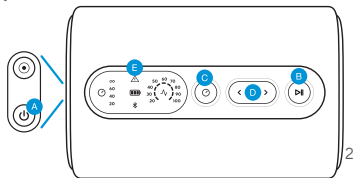2

## Getting started

1. To turn on the RecoveryAir Prime pump, press and hold the power button.
2. To get started, put your boots on first before attaching them to the pump. This allows you to get into the most comfortable position and avoid any disruption in the air flow connection.
3. Next, plug in the connector from the boots into the pump. Make sure you hear a “click” so that you know the RecoveryAir Prime is fully connected, is plugged in, and there is a proper seal.
4. Once connected and comfortable, you are ready to customize your treatment.

## Prepare your treatment

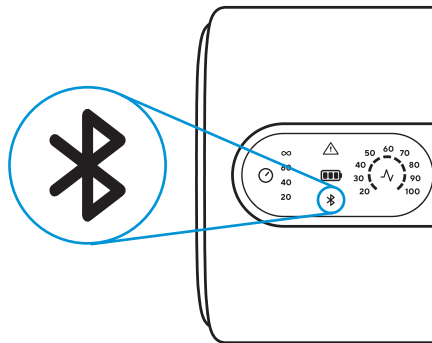
1. To turn on the RecoveryAir Prime pump, press and hold the power button.
2. Access two preset program options when connected to the Therabody app: Recover and Warm-up.
3. Once you have chosen your desired program, use the (<) and (>) buttons to adjust the pressure setting.
4. Then, toggle through the time options using the Time button to set your desired treatment time (20min, 40min, 60min, or continuous).
5. Once you have selected the treatment options that are right for you, press start, relax, and you'll be on your way to recovery in no time.

# Smart Features

The RecoveryAir Prime also connects to the Therabody app via Bluetooth, giving you control of your RecoveryAir Prime session from the convenience of your smartphone.

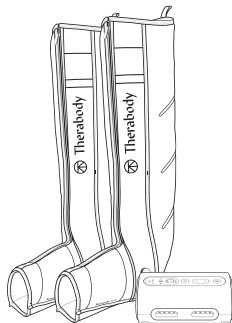
Not only can you easily access the Recover and Warm-up treatment options, but you can also adjust the pressure and time.

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